

MEBER⁺

MAX

POWERLIFT



POWER. SIMPLY POWER.

Hydraulically driven Power Stretcher
and Power Load system

www.meber.it


POWERLIFT

MAXIMUM LOAD

353Kg**778**lbs

**PUSH A BUTTON: LIFT. PUSH. DROP.
PUSH: LOAD. PUSH: DOWNLOAD.
THAT'S ALL, FOLKS.**

Max PowerLift System is the new combined stretcher/power load motorised solution. Max Powerlift was designed to simplify patient rescue and transport. Simplicity and immediacy are the core strengths of this device: every operation is characterised by ease of use, as simple as pressing a button, with the aid of a display that provides all info relevant to its functioning

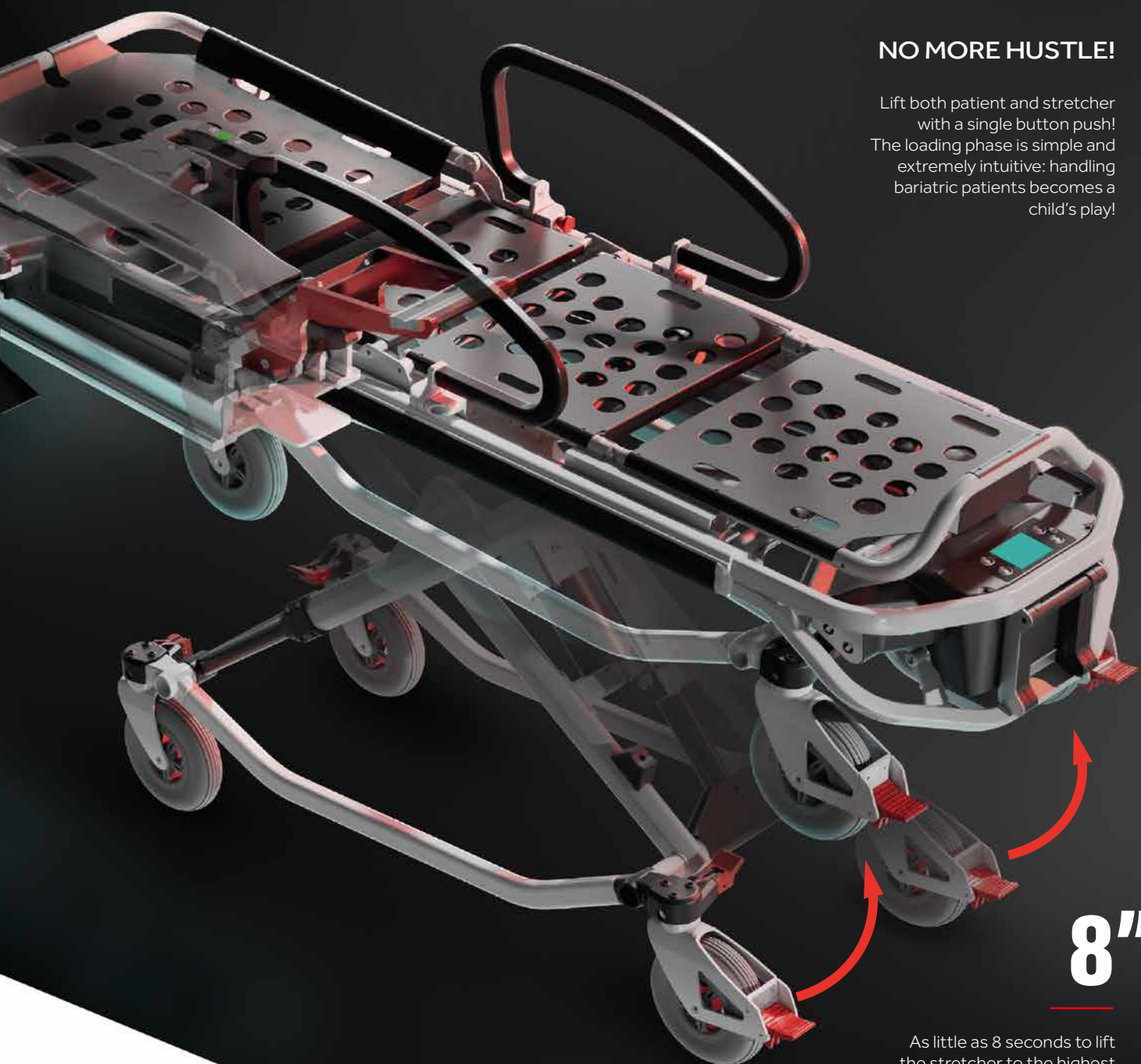

NO MORE HUSTLE!

Lift both patient and stretcher with a single button push! The loading phase is simple and extremely intuitive: handling bariatric patients becomes a child's play!

8"

As little as 8 seconds to lift the stretcher to the highest loading position!

MEBERT 

A super-powered-stretcher to make a difference at every despatch.

LED LIGHTS AND REFLECTIVE-, HIGH-VISIBILITY GRAPHICS

Max PowerLift is equipped with 4 LED lights, making it safe and easy to rescue patients in low-visibility conditions.



ALL HEIGHTS AT YOUR FINGERTIPS, WITH THE TOUCH OF A BUTTON. **AND MUCH MORE.**

Max PowerLift was designed to ensure safe and easy rescue operations. The hydraulically driven system allows to lift and drop the height of the stretcher all while being controlled by a battery-charged circuit.

FOWLER

Available Fowler position up to 30°.

TELESCOPIC FRAME

The telescopic head-side of the stretcher's frame can be retracted to add manoeuvrability in narrow spaces.

4 SWIVELLING CASTORS

Lock/Unlock systems for the 4 swivelling castors can be found on each individual wheel, with a self-aligning solution on the rear wheels.

BACKUP SYSTEM

The stretcher can be safely operated even in absence of its power supply




MAX POWERLIFT SYSTEM.

EXTRAORDINARILY SIMPLE. SUITABLE FOR ALL RESCUERS.

Simple. Extraordinarily simple... But how?

Max PowerLift factors in all state-of-the-art technology developed by Me.Ber., installed on a single device. To make it useful and usable, we resorted to a few expedients that simplified the usage of its system. Simplicity entails fewer training hours while optimizing efficiency. Max PowerLift adopts very intuitive commands on the sides of the stretcher and on the Power Load, with a single function per button.

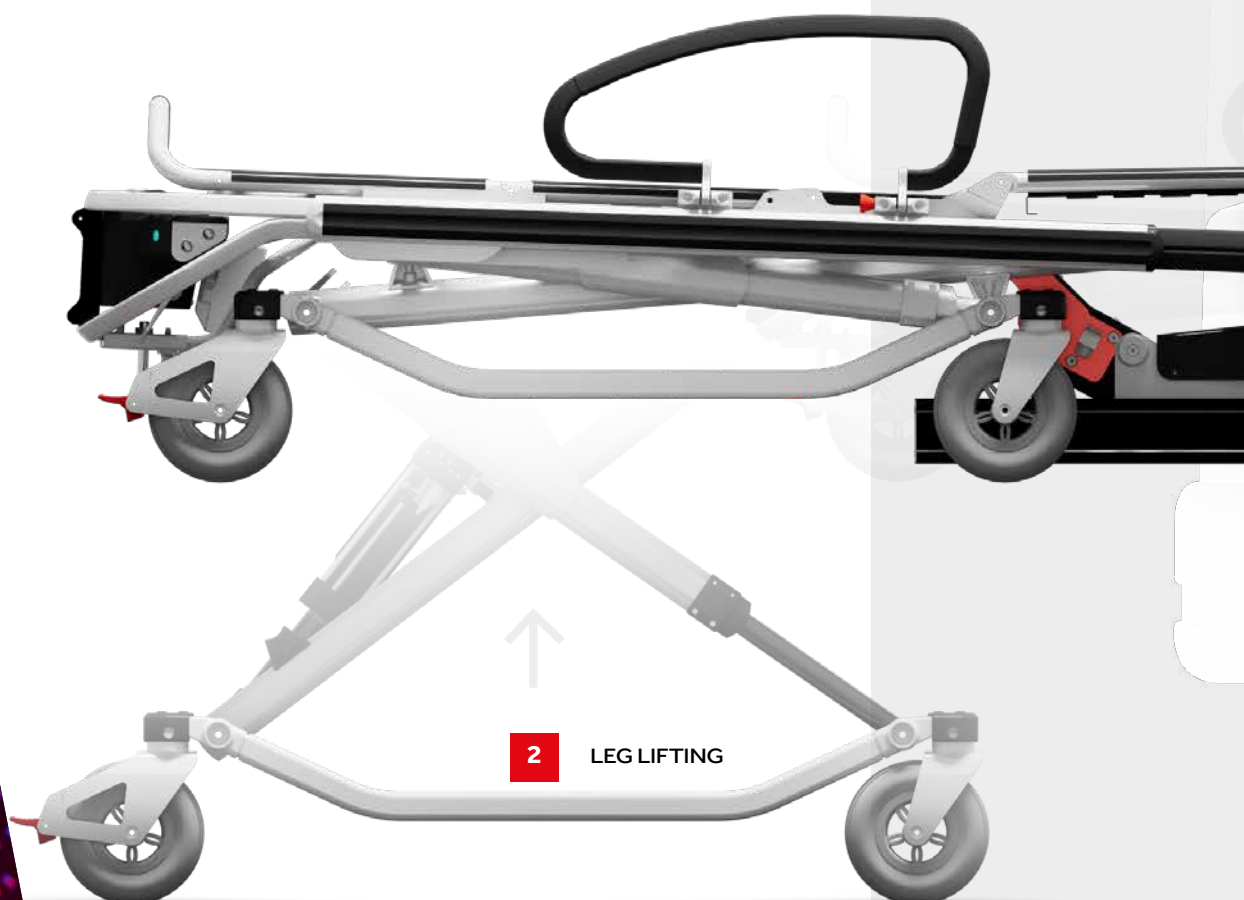

ADAPTABLE LOADING HEIGHT

The loading height can be customized, based on the vehicle.

ALWAYS AT THE RIGHT HEIGHT

You can adjust the Max PowerLift's height at will, continuously. Depending on the rescue scenario you may adjust the height to allow for ground rescue, then raise it up to the correct loading height, passing through countless intermediate positions.

3 LOADING

4 FASTENING

1 LOCKING

Hydraulically driven Power Stretcher and Power Load system



Max Powerlift Hydraulically driven Power Stretcher		
	X-frame ambulance cot	Loading System
Weight	75 kg	116 kg
Dimensions with extended trolley and extended head section	2150x595 mm	2100x695x240 mm
Material	Aluminum	Aluminum / Steel
Power supply	28V rechargeable battery	self recharging battery pack by ambulance power system
Max loading height	890 mm	
Backrest	reclining backrest from 0° to 85°	
Trendelenburg	0°÷15°	
Fowler	30°	
Loading capacity	350 kg in all heights	
From lowered position to maximum height	8"	
Wheels	4 ø200 swivel castors with directional lock	

All rights are reserved. Me.Ber. reserves the right to modify the content without prior notice. Me.Ber. s.r.l. Unipersonale cannot be held responsible for any damage caused by incorrect or incomplete information. Pictures shown are for illustration purpose only.

Accessories



Article Code
7002

ELECTRO-WELDED 4 PIECES
BLACK MATTRESS



Article Code
12128

ADJUSTABLE SYSTEM
FOR HEADREST ART.
12100/N



Article Code
12100/N

BLACK PILLOW FOR
HEADREST



Article Code
12092

BATTERY MILWAUKEE M28
RED 5,0 AMPERE.



Article Code
12090

BATTERY CHARGER MILWAUKEE
M28C FOR RECHARGEABLE
BATTERY M28 RED




Article Code
12120

LOWER - BLACK OBJECTS
HOLDER (SIZE: 92X43 CM)



Article Code
12122

UPPER - BLACK OBJECTS
HOLDER (SIZE: 39X52 CM)

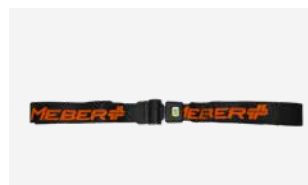

Article Code
Art. 608-N/MEB

BLACK 2 PIECES BELT WITH
METALLIC BUCKLE (L. 155 CM).



Article Code
612/N-2/MEB

ADJUSTABLE THORACIC
ABDOMINAL BELT BLACK OR
ORANGE/BLACK.

MEBER

Me.Ber. s.r.l. Unipersonale

Strada Langhirano, 270 43124 Fontanini (Parma) Italia

Tel. +39 0521 648770 | Fax +39 0521 648780 | e-mail: info@meber.it