



7097/4RG PROOF Mercury Lite Cinque 7097/4RG Proof

Certified aluminium 5 levels self-loading stretcher with 4 swivel wheels - high version



FEATURES

Mercury Cinque Lite high version. Mercury Cinque is derived from one of the most successful projects in the field of rescue, certainly the Me.Ber product most appreciated, the Mercury 4RG. Mercury Cinque Lite is intended to further extend the operational capabilities of the original model by introducing 5-level logic, five heights, as they are commonly called. Mercury Cinque Lite is much more than an evolution though! It is the only five-height with a folding front trolley (on the other hand, what do you need a stretcher to help you load a patient from his bed, if you can't get through the narrow spaces of a house?) also the only one that guarantees you the maximum capacity of 250 kg on all five levels. Furthermore it is even lighter than the stainless steel model.

FIVE LEVELS

LOAD LEVEL - Loading height

Loading on the vehicle, handling on the ground, final transfer to the emergency room or operating room, treatment or triage on the site.

HIGH LEVEL - High rescue

Suitable as a resting position. Ideal height for the placement of a spinalized patient.

INTERMEDIATE - High bed

Typical height of a hospital bed or height suitable for seating ambulatory patients. Ideal height to transfer the patient from stretcher to chair and vice versa.

LOW LEVEL - Bed or sofa

Typical height of a sofa or a bed. This height allows easy transfer of a collaborating patient and simplifies the transfer on the stretcher.

BASE LEVEL - Ground level rescue

At only 30 cm from the ground, it allows all conventional operations with a maximum load of 250 kg.

Heights: 100 - 90 - 81 - 68,5 - 39 cm

Mean features:

- - 250 kg Maximum load capacity on all heights
- Folding front trolley
- Ergonomics for back pain: the stretcher controls and general ergonomics are designed not to overload the rescuer, especially to relieve the load on the back.
- 200mm wheels to overcome any roughness. Greater comfort for the patient and great safety for rescuers who can move easily even in the presence of uneven ground and surfaces.
- 4 swivel wheels for safe manoeuvring at all times. Mercury Cinque has all four swivel wheels. A simple command allows the unlocking of the rotation of the front wheels for maximum manoeuvrability.


- Reflective decorations to improve the visibility of the stretcher in the roadside assistance context at night, reflective stickers were used in various points.

Compatible with fixing system art. 955 Proof.



Device certified according to the European harmonized safety standards UNI EN 1865 - UNI EN 1789.

MAIN FEATURES

10G Certification	
Material	Painted aluminium
Length (cm)	197
Length with folded trolley (cm)	167
Width (cm)	57
Folded height (cm)	31
Max loading height (cm)	71
Min loading height (cm)	69
Weight (kg)	42
Loading capacity (kg)	250

IMAGES



SEE ALSO



Mercury Lite Cinque 7095/4RG PROOF
Certified aluminium 5 levels self-loading stretcher with 4 swivel wheels



Mercury Cinque 7075/4RG PROOF
Certified 5 levels self-loading stretcher with 4 swivel wheels



Mercury Cinque 7077/4RG PROOF
Certified 5 levels self-loading stretcher with 4 swivel wheels - high version

ACCESSORIES



12080

"MPB" wide bariatric platform for self-loading stretchers Frog Mercury and Gecko complete of belts with reel and fixing for stretcher



12080-001

Monting kit for bariatric platform on self loading stretcher



12106

Mattress for AIR-COM stretchers - Anti-bedsores mattress with alternate segments



955 PROOF

EN 1789 fixing system for Proof roll-in stretchers



12154

"BEAR PODS" Bariatric containment system for self-loading stretchers



12018N

Telescopic foldable phleboclysis stand for self loading stretchers



939

Instrument support for stretchers



12128

Adjustable system for headrest of MeBer-X and Mercury stretchers



12100/N

Black electrowelded shaped cushion with MeBer logo for electrowelded mattresses



12070/BG

"COVER OX 2-3 BG"- professional cylinder holder for 2 and 3 lt oxygen bottles



12070/BR

COVER OX2-3 BR professional cylinder holder with belt and orange MeBer logo for 2 and 3 lt oxygen bottles



12072/BG

"COVER OX 5 BG"- professional cylinder holder for 5 lt oxygen bottles



12072/BR

"COVER OX 5 BR"- professional cylinder holder for 5 lt oxygen bottles



12084-001

Cover for bariatric wings set



12170

Black mattress for right bariatric wing for self-loading stretchers



12172

Black mattress for left bariatric wing for self-loading stretchers



7013/A
3 pcs. Skay fire-proof electrowelded
mattress



7014/A
4 pcs. Skay fire-proof electrowelded
mattress



938
Oximeb oxygen cylinder holder



12054
Oximeb Maxi - Oxygen cylinder
holder



707
Red 6-handle fireproof cover with
Me.Ber tape for self-loading
stretcher



608/A-MEB
2 pieces orange belt with metallic
hook and MeBer band (L. 155 cm)



608/N-MEB
2 pieces black belt with metallic
hook and MeBer band (L. 155 cm)




612/AN-2/MEB
Adjustable thoracic abdominal belt
orange/black with MeBer band and
4 points fastening for self loading
stretchers and sedan chairs



612/N-2/MEB
Adjustable thoracic abdominal belt
black with MeBer band with 4
points fastening for self loading
stretchers and sedan chairs



12025
Black multi usage object holder for
automatic loading stretcher

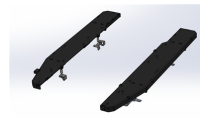

12200
Bariatric wings set for Mercury self-
loading stretchers complete with
mattresses, belts with reel and
assembly kit

STANDARD ACCESSORIES



12025
Black multi usage object holder for
automatic loading stretcher



608/N-MEB
2 pieces black belt with metallic
hook and MeBer band (L. 155 cm)



612/N-2/MEB
Adjustable thoracic abdominal belt
black with MeBer band with 4
points fastening for self loading
stretchers and sedan chairs



7013/N
3 pcs. Skay fire-proof electrowelded
mattress