



## 16150/V Ergon

---

Variable geometry comfort stretcher - green/blue



---

### FEATURES

[Visit Ergon website](#)

**ERGON** was designed and developed by rescue operators and experts of the medical emergency sector to offer outstanding performances, never achieved before. This stretcher ensures maximum comfort, especially in bed sore areas, thanks to a special material (the exclusive Me.Ber. **HC<sup>3</sup>** - High Comfort Composite Cell) consisting in a soft material moulded on a high resistance plastic polymer.

Besides patient's comfort, **ERGON** is designed to assist rescuers and operators thanks to unique solutions. The handles are upholstered in **HC<sup>3</sup>** and present a 20° inclination, thus allowing a better grip. Dedicated eyelets and supports were integrated for the fastening of straps and carabiners. When split, the distance between the two halves measures 25 mm: this ensures excellence support to the spinal column, in contrast with regular scoop stretchers which presents an 80 mm distance. The shape of the blade is also larger than usual and is easily placed between the patient and the ground. The employment of **HC<sup>3</sup>** material, besides conferring excellent structural stability, ensures great thermal insulation. This feature prevents hypothermia and allows for a longer stay of the patient on the Ergon stretcher, thanks to its design.



A technological revolution: the new jointed locking system FlexiLock allows to separate the stretcher in two halves and facilitates rescue operations even in extreme conditions. Besides fitting all standard rescue techniques, **ERGON** allows to easily perform rescue operations in extreme environments and when the patient is in critical conditions.



---

## MAIN FEATURES

---

Size	Adults
Usage	Rescue Immobilization Transport
Height (cm)	7,5
Width (cm)	43
Length (cm)	165
Max length (cm)	201
Folded length (cm)	120
Epaisseur pliée (cm)	9
Weight (kg)	10,5
Loading capacity (kg)	170
Trauma	
Alpine	
Material	HC3 - High Comfort Composite Cell Aluminium Aluminium alloy

---

## IMAGES

---



---

## SEE ALSO

---



**Ergon 16150/G**  
Variable geometry comfort  
stretcher - yellow/blue

---


---

## ACCESSORIES

---



**X-Lock 9150**  
Ergonomic head immobilizer for  
Ergon stretcher

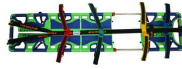

**X-FIX 9156**  
Restraint system with straps and  
carabiners



**959**  
Wall mount system for variable  
geometry stretcher Ergon



**16150-002**  
Storage bag for Ergon stretcher



**16150-001**  
Winch bag for Ergon



**584**  
Black bag with 3 orange straps in 2  
pcs and quick-release plastic  
buckle for Ergon stretcher



**572**  
Blue strap in 2 pieces with metallic  
buckle (l. 155 cm)



**574**  
Blue bag with 3 belts art. 572

---

## STANDARD ACCESSORIES

---



**584**  
Black bag with 3 orange straps in 2  
pcs and quick-release plastic  
buckle for Ergon stretcher

---