




9010 Galaxy

"GALAXY" Yellow spinal board

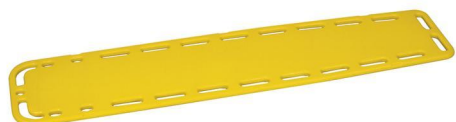

FEATURES

Spinal board made of plastic material at high strength. Designed paying particular attention to the difficulties met during every day's interventions; they have 4 holes for the quick and total fixing of the head immobilizer and two cavities where the board lays on the floor, when the base is blocked in the traditional way, that allow to avoid damages to the rip-off straps during the usage of the spinal board or the accommodation in the ambulance. It has 12 handles for the transport and is supplied with 3 belts with rapid unhooking buckle (Art. 605/MEB). X-ray translucent and available in 3 different colours: yellow, orange and red.

Device certified according to the European harmonized safety standards UNI EN 1865

MAIN FEATURES

Size	Adults
Length (cm)	183
Width (cm)	41
Thickness (cm)	5
Weight (kg)	5,5
Loading capacity (kg)	150
Trauma	
Usage	Immobilization
Material	ABS

SEE ALSO



Galaxy 9016
"GALAXY" Red spinal board



Galaxy 9012
"GALAXY" Orange spinal board

ACCESSORIES



Jelly 629

Orange standard head immobilizer
for spinal board



12144/N

Carrying case for spinal board



Spider 690

Belts immobilization system Spider
for spinal board

STANDARD ACCESSORIES



605/MEB

50 mm orange belt in one piece
with plastic buckle and ME.BER.
ribbon
