



## 683 Proxima Pro

Fixed padded sedan chair





### FEATURES

Made of light aluminium alloy, it differs from other "Proxima" versions for its light weight. Back rest and seat are padded with soft fireproof material covered by skay, granting high comfort and easiness to clean. Supplied with 4 wheels Ø 150 mm (2 rear rotating with stationing brake and 2 front fixed), 2 rear handles, 4 extractable levers (2 rear and 2 front, Ø 35 mm for an effective sliding) and foldable arms covered by an high comfort padded sheath. The wheels could be swapped by placing the directional brake on the back side and the stationing brake on the front side. Supplied with thoracic- abdominal belt. The UNI EN 1789 certified dedicated fixing system Art. 620 is available as accessory (see "Fixing systems").



Device certified according to the European harmonized standards UNI EN 1865 - UNI EN 1789


### MAIN FEATURES

10G Certification	
Wheel type	2 Rotating Wheels + 2 Small Fixed Wheels
Type	
Height (cm)	118
Width (cm)	57
Length (cm)	90
Max length (cm)	163
Weight (kg)	21
Loading capacity (kg)	170
Cardio	
Material	Painted aluminium

---

## IMAGES

---

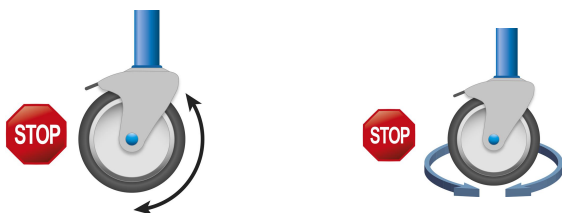

---

## ACCESSORIES

---



**988**  
Integrated control fastening system  
for Extra Ergolift FX and Proxima  
Pro chairs (optional)

---

## STANDARD ACCESSORIES

---



**610**  
Thoracic-abdominal belt with four  
fastening for padded chair

---